

# PROFESSIONAL URBAN FIRST AID



SAVE LIFE IN 4 MINUTES  
AED(Automated External Defibrillator)



SIMPLY TO USE



QUICK SHOCK



24/7/365 AUTOMATIC SELF-TEST



CLEAR AUDITORY AND VISUAL GUIDANCE



yuwell | PRIMEDIC  
Saves Life. Everywhere.

# AED



Primedic  
HeartSave Y | YA



**PRIMEDIC**  
Saves Life. Everywhere.



■ **Operation Modes**

Fully Automatic and Semi-Automatic models available

■ **Display Options**

Choose from models with no display, a 4.3" screen, or a 7" screen

■ **Multilingual Support**

Preloaded with 4 languages, selectable from a total of 32 options

■ **Energy Escalation Levels**

Adults: 200J → 300J → 360J  
 Children: 100J → 150J → 170J

• **Data Transfer**

Enables secure and reliable data retrieval via USB mobile storage, providing a convenient solution for users.

• **Design**

Features a fully sealed construction with no exposed interfaces, ensuring enhanced durability and protection.

• **Cold Environment Performance**

Operates reliably for a minimum of 20 minutes after transitioning from room temperature to -20°C.

• **Protection Standard**

Certified IP55 water- and dust-proof, making it suitable for use in various environments.

**30:2**

Support multiple compression modes of 30:2 and 15:2



Provide cycle sound guidance for CPR



Real-time feedback of CPR operation quality (optional)



Voice encouragement during CPR

## Steps for Using the AED



1

Open the Over



2

Attach Pads



3

Follow the Voice / Interactive Display



**Semi-Automatic:**

Click the Shock Button Required for Shock



**Fully Automatic:**

Automatically Shocks

## ■ Quick Operation



Turn on immediately by opening the cover, the operation interface is focused with visualization, less thinking time required.



# To Save Lives Everywhere

Cardiac arrest can happen to anyone, anywhere, anytime

Sudden Cardiac Death (SCD) is the abrupt cessation of the heart's pumping function, leading to an unexpected interruption of blood flow to the brain and vital organs. This results in loss of consciousness, respiratory arrest, and sudden death without treatment.

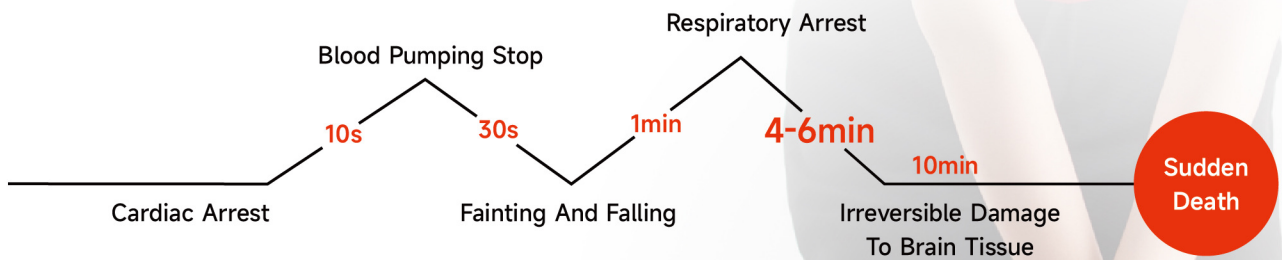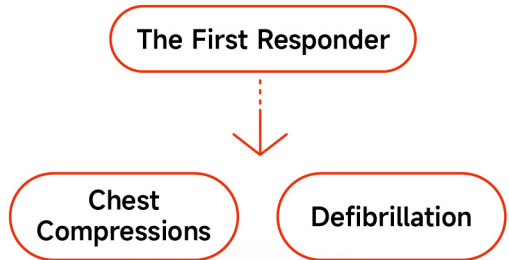

Saves Life,  
Everywhere

The PRIMEDIC Heartsave Y/YA series AED features advanced intelligent biphasic waveform defibrillation technology, ensuring a high first-attempt success rate with minimal myocardial damage, which automatically adjusts energy levels based on the patient's impedance, making it suitable for a wider group.

Brand	Primedic	Prompt Type	Voice / Interactive Display
Operation Mode	Fully & Semi Automatic	Warranty	5 years
Water/Dust Resistance (IP Rating)	IP55	Electrode Pads Life	5 years
Shock Energy	200-300-360(Adults) 100-150-170(Children)	Languages	Preloaded with 4 languages, selectable from a total of 32 options
Voice Prompt	Yes	LED Prompt	Choose from models with no display, a 4.3" screen, or a 7" screen
Battery Type	LiMnO <sub>2</sub> ,24V, 1.2Ah, non-rechargeable	Pediatric Capability	Yes
CPR Feedback	Optional (Depends on Model)	Pad Type	Disposable defibrillation electrodes
Dimensions and Weight:	22 x 9.7 x 29.6 cm (L x W x H)	Weight	Approx.2.5kg (with energy module and pads)

## PRIMEDIC From Germany

Focusing on defibrillation  
for more than 50 years

PRIMEDIC defibrillators are independently developed and designed by Metrax GmbH, headquartered in Rottweil, Germany. As the only specialized defibrillation equipment manufacturer with over 50 years of dedicated expertise, Metrax holds third-generation defibrillation patented technology and distributes its products across 128 countries.

The equipment is recognized for its safety and reliability, having obtained ISO 13485:2003 international quality system certification and CE certification.